


Qilive

Electric toothbrush Q.5194

- FR** Brosse à dents électrique
- ES** Cepillo de dientes eléctrico
- PT** Escova de dentes elétrica
- PL** Elektryczna szczoteczka do zębów
- HU** Elektromos fogkefe
- RO** Periuță de dinți electrică
- RU** Электрическая зубная щетка
- UA** Електрична зубна щітка

EN	User manual	P. 5
FR	Manuel d'utilisation	P. 15
ES	Manual de instrucciones	P. 25
PT	Manual de utilização	P. 35
PL	Instrukcja obsługi	S. 45
HU	Használati utasítás	55. O.
RO	Manual de instrucțiuni	P. 65
RU	Руководство пользователя	C. 75
UA	Довідник користувача	C. 86

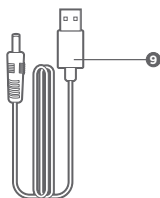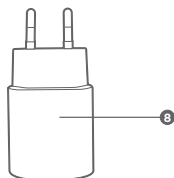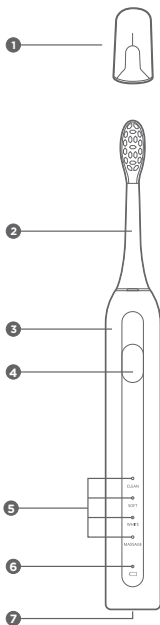EN

FR

ES

PT

PL

HU

RO

RU

UA

EN

FR

ES

PT

PL

HU

RO


RU

UA

1

2

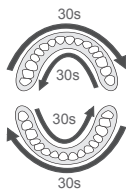
3


4

5

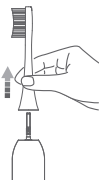6

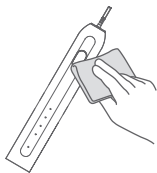CONTENT OVERVIEW :

EN

1. SAFETY INSTRUCTIONS	P. 5
2. TECHNICAL SPECIFICATIONS	P. 9
3. PRODUCT DESCRIPTION	P. 10
4. BEFORE FIRST USE	P. 11
5. CHARGING	P. 11
6. OPERATION	P. 11
7. CLEANING AND CARE	P. 13
8. TROUBLESHOOTING	P. 14

1. SAFETY INSTRUCTIONS

Before using the appliance, read this instruction manual carefully and keep it in a safe place for future reference.

DANGER – To reduce the risk of death or injury by electric shock:

- Do not place, store, or charge the appliance where it can fall or be pulled into a tub or sink. If the appliance falls into water, unplug it immediately first. Do not reach into the water.
- Do not immerse the appliance in water or any liquid.
- Keep the charger away from water.
- Make sure the appliance is completely dry before connecting it to the charger and power outlet.

WARNING – To reduce the risk of fire, electric shock, injury, death, or damage to the appliance:

- Risk of choking or suffocation! This product contains small parts. Keep the appliance and its cord out of reach of babies and children.
- This appliance is not intended to be used by children under the age of 8.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the

hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

- Close supervision is required when this product is used by or near children.
- This appliance is intended for personal and household use only.
- This appliance is not intended for medical or professional use on multiple patients in dental clinics, practice institutions, or similar settings.
- This appliance has been designed exclusively for brushing teeth.
- Do not use the appliance for any purpose other than its intended use as described in this manual.
- The appliance is designed for indoor use only.
- Do not leave the appliance unattended while it is charging. Unplug the appliance after charging and before cleaning it.
- To clean the appliance, follow the CLEANING AND CARE section.
- Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- The appliance is only to be used with the power supply unit provided with the appliance.
- This appliance contains batteries that are non-replaceable.
- **WARNING:** For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance.
-  Batteries must not be disposed of with normal household waste. Recycle them at the collection points provided for this purpose. Contact the local authorities or your dealer for recycling advice. Correct disposal of old batteries helps protect the environment and human health.

-  The marking of this device (symbol of a crossed-out dustbin) complies with European Directive 2012/19/EU relating to "used electrical and electronic equipment". This marking means that this equipment, at the end of its useful life, cannot be disposed of with other household waste. The user is required to deliver it to operators who collect used electrical and electronic equipment. Collection operators, including local collection points, shops and local authorities, create an appropriate system for collecting equipment. Appropriate collection of used electrical and electronic equipment helps avoid harmful impacts on human health and the environment, resulting from the presence of hazardous components, as well as from improper storage and treatment of this equipment.

CAUTION:

- Do not use the appliance while bathing or showering, as steam can enter and damage the appliance.
- If your toothpaste contains peroxide, baking soda, or bicarbonate (common in whitening toothpastes), clean the brush head and handle thoroughly to prevent the plastic from cracking.
- Keep the appliance and cord away from heated surfaces and any source of heat. Exposing the appliance to high temperatures can overheat the rechargeable battery and cause leakage, fire, or damage to the appliance.
- Never insert any metal or foreign object into the openings of the appliance, as this can cause a short circuit or fire.
- Do not short circuit, incinerate, crush, or damage batteries to avoid the risk of leakage, fire, explosion, and/or injury.
- Only use the appliance with parts and accessories recommended by the appliance manufacturer. The use of accessories or attachments not recommended by the appliance manufacturer may cause injury, fire, or damage to the appliance.
- Never charge the appliance with a damaged cord or plug.

- If the appliance is not working properly, is broken or damaged in any way, stop using it.
- This appliance has no user-serviceable parts. Do not attempt to dismantle, repair, or modify the appliance. This appliance has a built-in Lithium-ion battery and recyclable materials. Dispose of the appliance safely at your local recycling centre to help prevent soil, water and sea contamination.

MEDICAL WARNINGS

- If you have had oral or gum surgery in the last 2 months, consult your dentist before using this appliance.
- Consult your dentist if excessive bleeding occurs after using this appliance or if bleeding continues to occur after 1 week of use.
- If you have a pacemaker or any implanted medical device, consult your doctor or the manufacturer of the medical device before using this appliance.
- Stop using the appliance and consult your doctor or dentist if you experience any discomfort or pain.

2. TECHNICAL SPECIFICATIONS

Charger	
Input:	100-240 V~, 50/60 Hz
Output:	5 V  0.6 A
Output power	3.0 W
Manufacturer's name or trademark, commercial registration number and address	Shenzhen Tianyin Electronics Co., Ltd. 91440300788337003T 1st-6th Floor, Building 13, Tangtou Nangang third industrial zone, Tangtou community, Shiyan Town, Bao'an District, Shenzhen, Guangdong, P.R. China
Model identifier	TPA-147K050060VU01
Average active efficiency	73.32% min. (at 230V/50Hz)
Efficiency at low load (10 %)	/
No-load power consumption	0.060 W max. (at 230V/50Hz)

Appliance	
Charging input:	5 V  0.6 A
Battery:	3.7 V, 650 mAh Li-ion rechargeable battery, non-removable
Charging time:	≤ 4 hours
Operating time:	40 days (4 minutes or 2 twice a day in soft or massage modes) 20 days (4 minutes or 2 twice a day in standard or whiten modes)
Rated power:	1.6 W
Frequency range:	4 MHz
Operating/Storage temperature:	0 ~ 40 °C / -10 ~ 60°C
Ingress protection:	IPX7

3. PRODUCT DESCRIPTION

- | | |
|-----------------------------------|----------------------------|
| 1. Travel cap | 6. Battery indicator light |
| 2. Brush head | 7. Charging input |
| 3. Handle (main unit) | 8. Charging cable |
| 4. Power and Mode button | 9. Charger |
| 5. Mode indicator light (4 modes) | |

4. BEFORE FIRST USE

- Unpack the appliance.
- Check that the appliance has no visible damage and that all parts are complete.
- Fully charge the appliance.

5. CHARGING

Charge the appliance for 4 hours before first use, or when it has not been used for more than 3 months. Thereafter, a full charge will take approximately ≤ 4 hours.

1. Charge the appliance on a flat, dry, and stable surface.
2. Turn off the appliance.
3. Connect the charging cable to the appliance and plug the charger into a suitable power outlet. The battery indicator light will flash red when the battery of the appliance is charging and will turn off when fully charged. Do not leave the appliance unattended while it is charging.
4. Once charging is complete, unplug the charger and appliance from the power outlet.
 - Recharge the appliance when the battery indicator flashes red.

6. OPERATION

CAUTION: Do not brush excessively or apply too much pressure while brushing as this can cause gum recession, tooth sensitivity, or enamel erosion. Do not use the appliance for more than the recommended time.

It is recommended to use the appliance twice a day, for 2 minutes each time.

1. Insert the brush head into the handle with the bristles pointing forward ①.

NOTE: A small gap between the brush head and the handle is normal for operation.

2. Wet the brush head and apply a small amount of toothpaste.
3. Press the Power button to turn on the toothbrush.
4. To switch the brushing mode, press the Power button repeatedly. The indicator light of the mode selected will turn on. Four brushing modes are available:
 - Clean mode: Standard mode for optimal cleaning.
 - Soft mode: Gentle cleaning for sensitive teeth and gums.
 - Whiten mode: Deep cleaning to polish and whiten teeth.
 - Massage mode: Gentle cleaning with vibrating pulses to stimulate and massage the gums.

NOTE: The toothbrush will pause briefly every 30-second interval to remind you to switch to another section and will automatically turn off after the recommended 2 minutes. The timer will restart if the mode is switched during brushing.

5. Divide your mouth into 4 sections and brush each section for 30 seconds ②. Start with the outer surfaces of the upper teeth, then continue with the inner and biting surfaces. Repeat on the lower teeth.
6. Place the brush head at a 45° angle toward the gumline ③. Apply light pressure and move the brush head along the gumline and surfaces of the teeth, guiding the bristles slowly from tooth to tooth with short back and forth strokes.
7. The appliance will automatically turn off and remember the mode used.
8. After the appliance is turned on, press the Power button immediately (within 3 seconds) to switch modes. When the appliance is turned on for more than 3 seconds, press the Power button once to turn it off.

7. CLEANING AND CARE

EN

WARNING: Before cleaning, always turn off and unplug the appliance from the power outlet.

WARNING: To protect against fire, electric shock, or damage to the appliance, never immerse the handle (main unit), charging cable or charger in water or any liquid. Keep these electrical parts dry.

1. After each use, rinse the brush head thoroughly under the tap with cold or lukewarm water ④. Remove the brush head from the handle to rinse them separately ⑤.
2. Dry the handle with a clean, dry towel ⑥. Remove excess water from the brush head. Allow to air dry before reassembling.
3. Make sure the handle is completely dry before charging.
4. Store the appliance in a cool and ventilated place. Keep away from direct sunlight.
5. Replace your brush head approximately every 3 months.
6. Never use solvents, chemical or abrasive cleaning agents, wire brushes, sharp objects or scouring pads to clean the appliance or accessories.

8. TROUBLESHOOTING

Problem	Possible cause	Possible solution
The appliance does turn on.	■ The battery is discharged. ■ The appliance has not been used for more than 3 months.	■ Fully charge the battery for 4 hours. ■ Avoid deep discharge to prolong battery life.
Weak or reduced vibration.	The toothbrush is pressed too hard against the teeth.	Place the toothbrush gently against the teeth.
The toothbrush operates only for a short period of time even immediately after being charged.	The battery has reached the end of its useful life.	Do not dispose the appliance with normal household waste. Recycle the appliance in your nearest recycling centre.

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SAS OIA, 200, rue de la Recherche, Le Colibri